	文件类型	技术文件	文件版本	V3.3.0
	项目代号	/	页 码	第0页 共 1页
	文件编号	GL DM-56-17-en		

文件名称：  
**CONTROL DM-RET RET0523 2023-07-10**

适用范围：1116 项目

编制：	审核：	批准：
日期：	日期：	日期：

### 修 订 记 录

版本	TCN/ECR/PCN 编号	修订内容概要	修改人	批准人	批准日期
V3.3.0	/	新建	黄江妮	郑剑通	2023.05.19


发至：  
总经办 质量部 财务部 营销部 客服 计划 采购 工程 仓库  
设备研发 试剂研发 物料质量 设备生产 试剂生产 试剂质检 设备质检


由试剂研发

# CONTROL DM-RET

## DM-RET

### HEMATOLOGY CONTROL


**CONTROL**
 2023-04-12

 2023-07-10


Applicable Instruments	Parameter	Unit	LOT RET0523L	LOT RET0523N	LOT RET0523H
DYMIND	RET#	$\times 10^{12}/L$	0.0368 $\pm 0.0185$	0.2243 $\pm 0.1316$	0.2671 $\pm 0.1200$
DH-600	RET%	%	0.65 $\pm 0.33$	5.02 $\pm 2.01$	11.41 $\pm 4.50$
DH-602	LFR	%	89.8 $\pm 27.0$	89.7 $\pm 27.0$	89.4 $\pm 26.9$
DH-605	MFR	%	8.1 $\pm 12.2$	8.2 $\pm 12.4$	8.6 $\pm 13.0$
DH-610	HFR	%	2.1 $\pm 3.3$	2.0 $\pm 3.1$	1.9 $\pm 3.0$
DH-612	IRF	%	10.2 $\pm 15.4$	10.3 $\pm 15.4$	10.6 $\pm 15.9$
DH-615 (Technical File Version A1.6 or higher)	RHE	pg	16.3 $\pm 10.0$	16.2 $\pm 6.5$	16.7 $\pm 6.7$

**【NOTE】**


1. The controls should be stored in refrigerator (2°C~8°C). After opening, it will keep stable for 14 days when it is stored airtight at 2°C~8°C.
2. Before mixing and running the control after take it out from the refrigerator, please keep it at least 15 minutes until reaching room temperature(15°C~30°C).
3. Controls must be well mixed before using. Please mix gently, to avoid cells rupture and/or generating bubbles.
4. After using, put the controls back into the refrigerator to prevent contamination and evaporation.

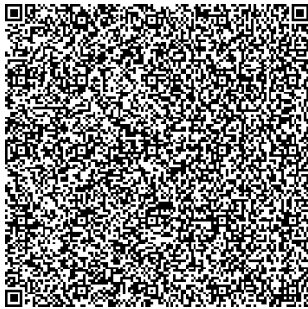

RET0523L

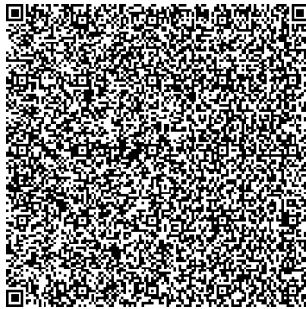

RET0523N

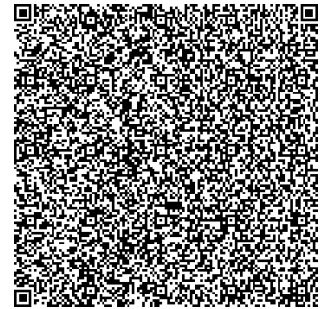

RET0523H