	文件类型	技术文件	文件版本	V3.6.0
	项目代号	/	页 码	第0页 共 2页
	文件编号	GL DM-56-17-en		

文件名称:

**CONTROL DM-RET RET1123 2024-01-10**

适用范围: 1116 项目、7003 项目

编制: 王雨婷

审核: 刘金盈

批准: 郑剑通

日期: 2023.10.31

日期: 2023.10.31

日期: 2023.10.31

修 订 记 录

版本	TCN/ECR/PCN 编号	修订内容概要	修改人	批准人	批准日期
V3.6.0	/	新建	王雨婷	郑剑通	2023.10.31

发至:

总经办 质量部 财务部 营销部 客服 计划 采购 工程 仓库

设备研发 试剂研发 物料质量 设备生产 试剂生产 试剂质检 设备质检

由试剂研发

# DM-RET

## HEMATOLOGY CONTROL

**CONTROL**


📅 2023-10-13

🕒 2024-01-10


Applicable Instruments	Parameter	Unit	LOT RET1123L	LOT RET1123N	LOT RET1123H
DYMIND	RET#	×10 <sup>12</sup> /L	0.0428 ±0.0214	0.2329 ±0.1397	0.2813 ±0.1406
DH-600	RET%	%	0.76 ±0.38	5.04 ±2.02	12.02 ±4.81
DH-602					
DH-605	LFR	%	90.0 ±27.0	89.5 ±26.8	89.9 ±27.0
DH-610					
DH-612	MFR	%	7.7 ±11.6	8.6 ±12.9	8.2 ±12.3
DH-615					
(Technical File Version A2.0 or higher)	HFR	%	2.2 ±3.4	2.0 ±3.0	2.0 ±2.9
	IRF	%	10.0 ±14.9	10.6 ±15.8	10.1 ±15.2
	RHE	pg	16.1 ±6.4	14.5 ±5.8	12.8 ±5.1

**【NOTE】**


1. The controls should be stored in refrigerator (2°C~8°C). After opening, it will keep stable for 14 days when it is stored airtight at 2°C~8°C.
2. Before mixing and running the control after take it out from the refrigerator, please keep it at least 15 minutes until reaching room temperature(15°C~30°C).
3. Controls must be well mixed before using. Please mix gently, to avoid cells rupture and/or generating bubbles.
4. After using, put the controls back into the refrigerator to prevent contamination and evaporation.

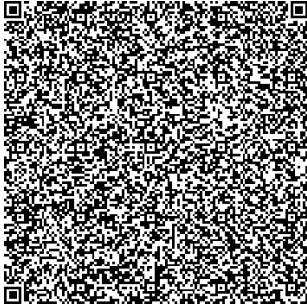

RET1123L

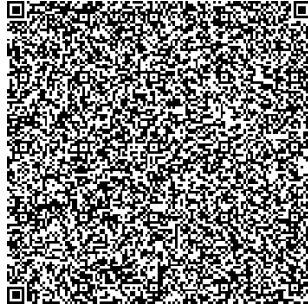

RET1123N

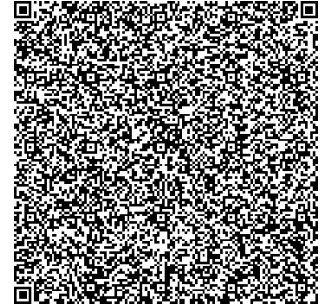

RET1123H

# DM-RET

## HEMATOLOGY CONTROL

**CONTROL**

📅 2023-10-13

🕒 2024-01-10

Applicable Instruments	Parameter	Unit	LOT RET1123L	LOT RET1123N	LOT RET1123H
DYMIND  DM60 Vet DM61 Vet DM62 Vet DM63 Vet (Technical File Version A1.1 or higher)	RET#	×10 <sup>12</sup> /L	0.0277 ±0.0222	0.1740 ±0.0870	0.1432 ±0.0716
	RET%	%	0.49 ±0.39	3.71 ±1.86	5.92 ±2.96

**【NOTE】**

1. The controls should be stored in refrigerator (2°C~8°C). After opening, it will keep stable for 14 days when it is stored airtight at 2°C~8°C.
2. Before mixing and running the control after take it out from the refrigerator, please keep it at least 15 minutes until reaching room temperature(15°C~30°C).
3. Controls must be well mixed before using. Please mix gently, to avoid cells rupture and/or generating bubbles.
4. After using, put the controls back into the refrigerator to prevent contamination and evaporation.



RET1123L



RET1123N



RET1123H