	文件类型	技术文件	文件版本	V4.2.0
	项目代号	/	页 码	第0页 共 2页
	文件编号	GL DM-56-17-en		

文件名称:

CONTROL DM-RET RET0324 2024-05-10

适用范围: 1116 项目、7003 项目

编制: 刘嘉 日期: 2024.02.29	审核: 刘金盈 日期: 2024.02.29	批准: 郑剑通 日期: 2024.02.29
--------------------------	---------------------------	---------------------------

修 订 记 录

版本	TCN/ECR/PCN 编号	修订内容概要	修改人	批准人	批准日期
V4.2.0	/	新建	刘嘉	郑剑通	2024.02.29

发至: 总经办 质量部 财务部 营销部 客服 计划 采购 工程 仓库
设备研发 试剂研发 物料质量 设备生产 试剂生产 试剂质检 设备质检

由试剂研发

DM-RET

HEMATOLOGY CONTROL

CONTROL


📅 2024-02-11

🕒 2024-05-10


Applicable Instruments	Parameter	Unit	LOT RET0324L	LOT RET0324N	LOT RET0324H
DYMIND	RET#	×10 ¹² /L	0.0430 ±0.0215	0.2962 ±0.1316	0.2497 ±0.1200
DH-600	RET%	%	0.75 ±0.37	6.53 ±2.61	10.90 ±4.36
DH-602	LFR	%	89.6 ±27.0	89.6 ±27.0	89.7 ±27.0
DH-605	MFR	%	8.2 ±12.0	8.2 ±12.0	8.5 ±13.0
DH-610	HFR	%	2.2 ±3.0	2.2 ±3.0	1.9 ±3.0
DH-612	IRF	%	10.4 ±16.0	10.4 ±16.0	10.3 ±16.0
DH-615	RHE	pg	16.8 ±7.0	17.1 ±7.0	17.3 ±7.0
(Technical File Version A2.0 or higher)					

【NOTE】


1. The controls should be stored in refrigerator (2°C~8°C). After opening, it will keep stable for 14 days when it is stored airtight at 2°C~8°C.
2. Before mixing and running the control after take it out from the refrigerator, please keep it at least 15 minutes until reaching room temperature(15°C~30°C).
3. Controls must be well mixed before using. Please mix gently, to avoid cells rupture and/or generating bubbles.
4. After using, put the controls back into the refrigerator to prevent contamination and evaporation.

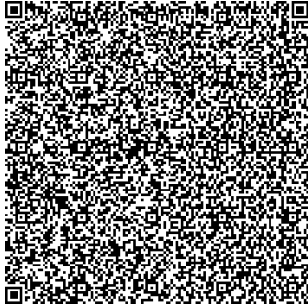

RET0324L

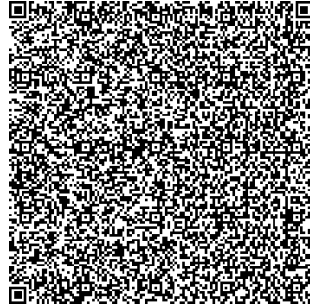

RET0324N

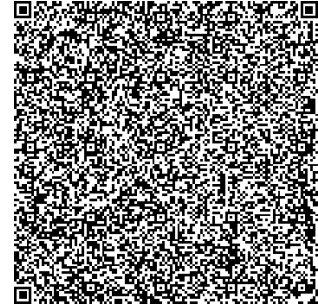

RET0324H

DM-RET

HEMATOLOGY CONTROL

CONTROL

📅 2024-02-11

🕒 2024-05-10

Applicable Instruments	Parameter	Unit	LOT RET0324L	LOT RET0324N	LOT RET0324H
DYMIND DM60 Vet DM61 Vet DM62 Vet DM63 Vet (Technical File Version A1.4 or higher)	RET#	×10 ¹² /L	0.0203 ±0.0162	0.2540 ±0.1270	0.1595 ±0.0798
	RET%	%	0.35 ±0.28	5.50 ±2.75	6.69 ±3.35

【NOTE】

1. The controls should be stored in refrigerator (2°C~8°C). After opening, it will keep stable for 14 days when it is stored airtight at 2°C~8°C.
2. Before mixing and running the control after take it out from the refrigerator, please keep it at least 15 minutes until reaching room temperature(15°C~30°C).
3. Controls must be well mixed before using. Please mix gently, to avoid cells rupture and/or generating bubbles.
4. After using, put the controls back into the refrigerator to prevent contamination and evaporation.



RET0324L



RET0324N



RET0324H