

# DH76

5-Part Auto Hematology Analyzer

<b>Principles</b>	Electrical Impedance method for WBC/BAS count and RBC/PLT tests Cyanide free colorimetry for HGB test Laser-based flow cytometry for WBC differentiation	<b>Sample volume</b>	20 $\mu$ L
<b>Parameters</b>	25 reportable parameters: WBC, Neu#, Lym#, Mon#, Eos#, Bas#, Neu%, Lym%, Mon%, Eos%, Bas%, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-SD, RDW-CV, PLT, MPV, PDW, PCT, P-LCR, P-LCC 4 research parameters: ALY#, ALY%, LIC#, LIC% 3 histograms for WBC, RBC and PLT One 3D scattergram and three 2D scattergrams for WBC differentiation	<b>Linearity range</b>	WBC: 0.00-300 $\times 10^9$ /L RBC: 0.00-8.5 $\times 10^{12}$ /L HGB: 0-250g/L PLT: 0-3000 $\times 10^9$ /L HCT: 0~67%
		<b>Repeatability</b>	WBC $\leq$ 2.0%(4.0-15.0 $\times 10^9$ /L) RBC $\leq$ 1.5%(3.5-6.0 $\times 10^{12}$ /L) HGB $\leq$ 1.5%(110-180g/L) MCV $\leq$ 1%(70-120fL) PLT $\leq$ 6.0%(100-149 $\times 10^9$ /L) $\leq$ 4.0%(150-500 $\times 10^9$ /L)
		<b>Communication</b>	LAN port supports HL7 protocol Support bi-directional LIS
<b>Printout</b>	External printer, compatible with multiply laser / inkjet printers, compatible with various formats and user-defined formats	<b>Sample mode</b>	Whole blood, capillary whole blood and pre-diluted modes
<b>Throughput</b>	Up to 90 tests per hour	<b>Power requirement</b>	100V-240V, 50/60Hz, $\leq$ 200VA
<b>Storage</b>	Up to 150,000 records	<b>Dimension</b>	650mm(W)x550mm(H)x610mm(D)
		<b>Net weight</b>	57.5kg

## SHENZHEN DYMIND BIOTECHNOLOGY CO.,LTD.

10th Floor, Building B, High-tech Park, Guangqiao Road, Tianliao Community, Yutang Street, Guangming District, Shenzhen 518107, P. R. China

+86-755-26008015-8123

Intl@dymind.com

www.dymind.com



Declaration: Shenzhen Dymind Biotechnology Co., Ltd reserves the right to change the product of specifications and appearance at any time. For the information of this manual, Shenzhen Dymind Biotechnology Co., Ltd reserves the right to the interpretation and the decision.  
P/N: EN-DH76[7.0]

 **DYMIND** | *the simpler, the better*

## Dedication to Efficiency



# DH76

5-Part Auto Hematology Analyzer

## Fully Automatic Intelligent Technology

Visual self-tests (working status self-test and reagent volume self-detection)

Automatic test re-run, with re-run rules inside

One-button error removal and automatic standby mode function

Intelligent WBC/RBC/PLT floating threshold technology




## Excellent Data Management System

A wide range of pathological alarms and up to 20 user-defined alarm ranges


Customized templates print and help to customize the final report according to requirements

Bi-directional LIS connection

QC and calibration management

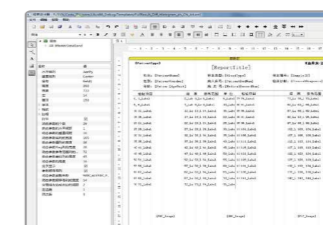

Especially for blood sample untreated over a long time

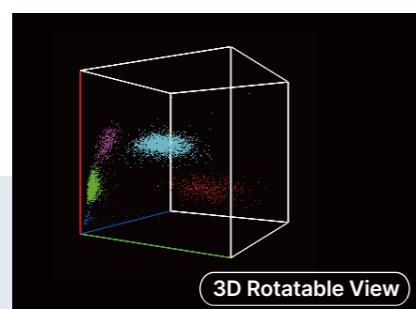


3D scattergram for accurate differentiation and pathological sample flag



Loading up to 60 samples simultaneously and supporting continuous samples loading



Up to 90 tests per hour, auto sample mixing

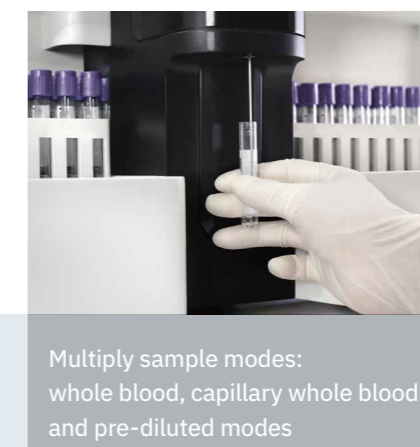

Multiply sample modes: whole blood, capillary whole blood and pre-diluted modes



Closed and open vial modes, emergency mode available